




ENTERAL NUTRITION  
Syringe pump



Syringe pump  
for enteral nutrition

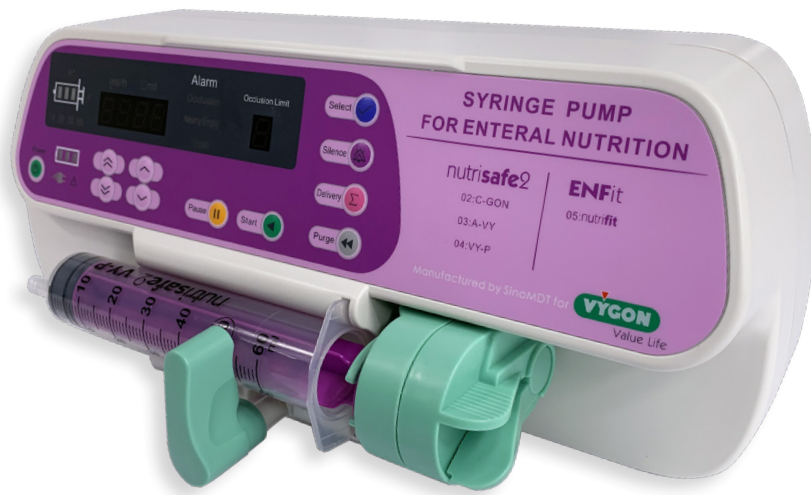

Value Life

## ● Accuracy in enteral feeding

Designed to meet all patients needs in stationary medical care, this easy-to-use high-quality syringe pump is dedicated to accurate enteral feeding.



The following syringes are already calibrated on the device. Others channels are available to calibrate new syringe brands (contact your biomedical department).

Syringe vol.	Nutrisafe 2			Nutrifit
	Channel 2 / CGON	Channel 3 / AVY	Channel 4 / VYP	Channel 5 / Nutrifit
10ml	1015.102	1015.103	1015.104	1015.107
20ml	1015.202	1015.213	1015.204	1015.207
30 / 35ml	1015.352	-	-	-
50ml	1015.602	1015.603	1015.604	1015.607

## Pump functions

- Automatic syringe size identification
- Optic and sound alarms
- Possibility to adjust sound level
- Up to 2000 history records
- Delivered volume inquiry
- Occlusion detection pressure setting
- Alarms: occlusion / prior administration end / administration end / battery low / out of battery / mains failure / equipment error / flowrate overlimit / administration completion / no operation within 2 minutes


## Advantages



For all brands of Nutrisafe 2 and Nutrifit syringes



6 hours battery autonomy

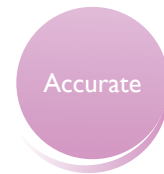

Accuracy of flowrate including syringe  $\pm 2\%$

## Technical features

Dimension (LxWxH mm)	306mm x 140mm x 146mm	
Weight	2.25kg (with clamp base)	
Accuracy of flow rate (Nominal value)	Accuracy including syringe: within $\pm 1\%$ Mechanical accuracy: within $\pm 1\%$	
Purge rate	1500 mL/h 900 mL/h 600 mL/h 400 mL/h	50 mL syringe 30/35 mL syringe 20 mL syringe 10 mL syringe
Flow rate	0.1 to 1500 mL/h (specified value 0.1 to 999.9 mL/h, 0.1 mL/h step) 0.1 to 900 mL/h (0.1 mL/h step) 0.1 to 600 mL/h (0.1 mL/h step) 0.1 to 400 mL/h (0.1 mL/h step)	50 mL syringe 30/35 mL syringe 20 mL syringe 10 mL syringe

## Product configuration

Product code: **SN-50C6**

- Syringe pump for enteral nutrition
- Power cords (Europe (plug C) / UK (plug G) / US (plug A))
- User manual and sticker
- Packing list
- Certificate of compliance
- Warranty card

## Contact

VYGON

8 rue de Paris 95440 Ecoeuven – France

Export customer service

Tel : + 33 (0)1 39 92 64 15

E-mail : [export@vygon.com](mailto:export@vygon.com)

Fax : + 33 (0)1 34 29 19 34

Service commercial France

Tel : + 33 (0) 1 39 92 63 81

E-mail : [savfrance@vygon.com](mailto:savfrance@vygon.com)

Fax : + 33 (0) 1 39 90 29 37

## OBSTETRICS NEONATOLOGY ENTERAL

**For further information, please contact: [marketingbenelux@vygon.com](mailto:marketingbenelux@vygon.com)**

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature.


---

SA Vygon NV • Haachtsesteenweg 1650 • Chaussée de Haecht 1650 • 1130 BRUSSELS • BELGIUM

Tel.: +32 (0)2 706 09 50 • [www.vygon.be](http://www.vygon.be)

Vygon Nederland BV • Kerkhofstraat 21 • 5554 HG VALKENSWAARD • THE NETHERLANDS

Tel.: +31 (0)40 208 93 80 • [www.vygon.nl](http://www.vygon.nl)

 Vygon Benelux

VY19163



[www.vygon.be](http://www.vygon.be)

[www.vygon.nl](http://www.vygon.nl)